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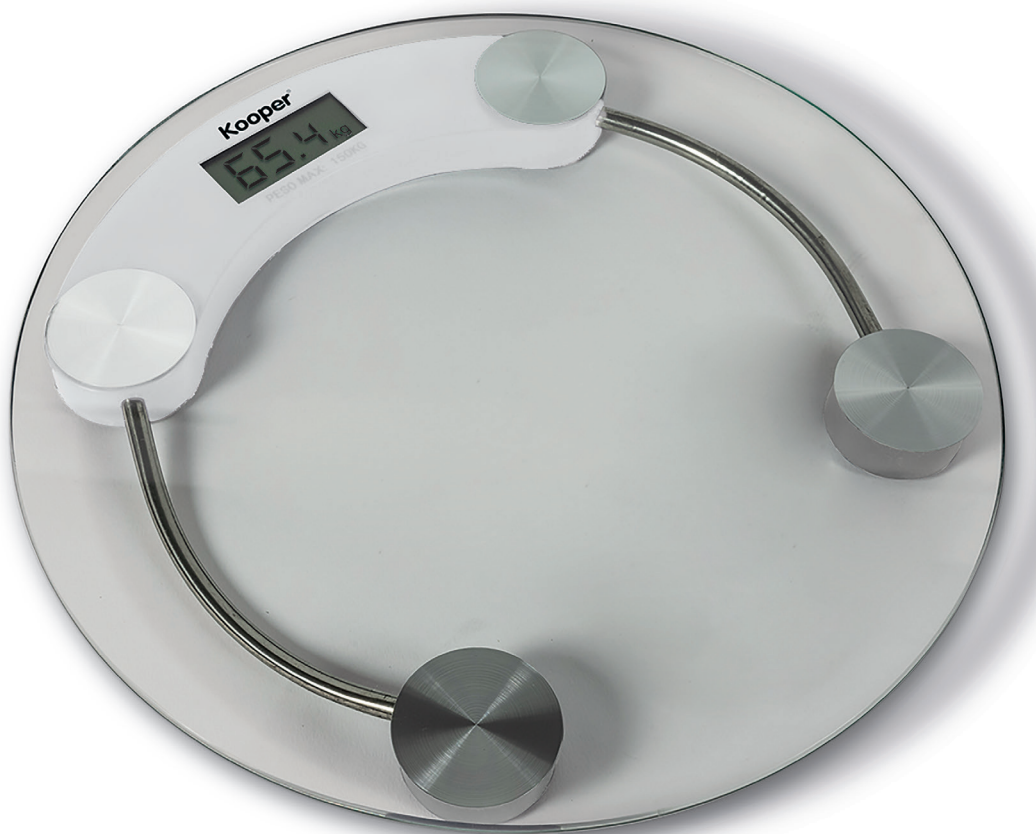
GRAZIE!

DISPOSAL OF THE ELECTRIC OR ELECTRONIC DEVICE AT THE END OF ITS LIFE

**INFORMATION FOR USERS ABOUT THE CORRECT DISPOSAL
OF THE PRODUCT ACCORDING TO THE EUROPEAN DIRECTIVE 2012/19/EU
AND TO THE LEGISLATIVE DECREE N. 49/2014.**

DEAR CUSTOMER,
IT IS IMPORTANT TO REMIND YOU THAT, AS A GOOD CITIZEN, YOU CANNOT DISPOSE OF THIS PRODUCT AS MUNICIPAL WASTE; YOU HAVE TO DO A SEPARATE COLLECTION. THE SYMBOL LABELLED ON THE APPLIANCE INDICATES THAT THE RUBBISH IS SUBJECT TO "SEPARATE COLLECTION". ELECTRICAL AND ELECTRONIC APPLIANCE CANNOT BE CONSIDERED AS DOMESTIC WASTE. THEREFORE, AT THE END OF ITS USEFUL LIFE, THIS PRODUCT SHOULD NOT BE DISPOSED WITH MUNICIPAL WASTES. THE USER MUST THEREFORE ASSIGN OR (HAVE COLLECTED) THE RUBBISH TO A TREATMENT FACILITY UNDER LOCAL ADMINISTRATION, OR HAND IT OVER TO THE RESELLER IN EXCHANGE FOR AN EQUIVALENT NEW PRODUCT.

THE SEPARATE COLLECTION OF THE RUBBISH AND THE SUCCESSIVE TREATMENT, RECYCLING AND DISPOSAL OPERATIONS PROMOTE PRODUCTION OF APPLIANCES MADE WITH RECYCLED MATERIALS AND REDUCE NEGATIVE EFFECTS ON HEALTH AND THE ENVIRONMENT CAUSED BY IMPROPER TREATMENT OF RUBBISH.

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ONLY FOR HOME USE. ATTENTION! READ THE INSTRUCTIONS GIVEN IN THIS HANDBOOK VERY CAREFULLY. IT SHOULD BE PUT ASIDE FOR FURTHER CONSULTATION DURING THE ENTIRE LIFECYCLE OF THE PRODUCT SINCE IT GIVES IMPORTANT INDICATIONS ON THE PRODUCT'S FUNCTIONS AND SAFETY DURING INSTALLATION, USE AND MAINTENANCE.

THIS ELECTRONIC PERSONAL SCALE HAS BEEN DESIGNED TO MONITOR BODY WEIGHT SPECIFICATION

Equiped with a high precision strain sensor system;
Graduation:100g/0.2LB;
Mini effective:0.5kg;
Automatic zero resetting;
Automatic switch off;
Unit:kg/Lb;
symbol Lo bat low-voltage detection, low voltage;
Hint of out:o-ld
Hint of wrong:---;
Sensor damaged or circuit breaker display:LCER;
Key press: unit conversion switch;program to load switch

How to weight

1. Remove the insulation film, put the scale on a hard & flat surface (avoid carpet of soft surface),wait. 2 seconds before weighing.
2. In order to change the unit of weight measurement (Kg or Lb), open the battery compartment and select the unit of interest using the switch located there
3. Set up unit prior to use, unit can net be converted when weight.
4. Stand on the scale in balance, you can read the result automatically when the weig ht is set after 3 seconds and automatically switch off after 8 seconds. .

Note: Before putting the scale in the flat & hard surface, the scale is automatically switch off, you must weight automatically first..

ADVICE FOR USE AND CARE

- When unused the scale should be left unloaded to avoid aver assumption of power and fai lure of work.
- Try to stand in the middle piace to avoid slip.
- Stand on the scale in good balance to avoid failure of lock and incorrect illustration.
- No heavy knock in the scale.
- Keep the scale dry. Do net use it in areas with water stains. Do net rinse with water.
- Handle with care, when using.

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BATTERY WARNINGS

- Remove the battery/battery pack/accumulator from the units if they are not to be used for long periods
- Do not recharge non-rechargeable batteries
- Rechargeable batteries/battery packs/accumulators must be removed from the product before charging them with an external charger.
- Charging must be supervised by an adult
- Insert batteries/battery packs/accumulators according to the correct polarity
- Installation and replacement of batteries/battery packs/accumulators must be carried out by an adult
- Do not immerse the item or batteries/battery packs/accumulators in water
- Do not expose batteries/battery packs/accumulators to heat or fire
- If liquid spilled from batteries/battery packs/accumulators comes into contact with clothing or skin, immediately wash the part with clean water and seek medical advice
- Stop using batteries/battery packs/accumulators immediately if you notice any anomalies
- Avoid contact of batteries/battery packs/accumulators with metal objects, such as keys or coins
- Batteries/battery packs/accumulators may explode or release toxic substances. Risk of fire or burns. Do not open, crush, modify, disassemble, pierce, heat or burn
- A battery pack is to be understood as a single unit which must not be torn, opened or, in any case, tampered with
- Do not put batteries/battery packs/rechargeable batteries in your mouth. If swallowed, promptly contact your doctor or local poison control centre
- Do not cause short circuits, to avoid burns or fire
- If batteries/battery packs/accumulators swell, discontinue use and dispose of devices in accordance with regulations
- Be very careful when using batteries/batteries/accumulators in series. Never mix batteries/batteries/accumulators with different characteristics in terms of brand, capacity, voltage. Do not mix partially discharged batteries/batteries/accumulators with charged batteries/batteries/accumulators.
- Ensure that batteries/batteries/accumulators are of the correct type and voltage for the appliance.

DANGER: DAMAGE AND/OR LEAKAGE!

FAILURE TO OBSERVE THE ABOVE MENTIONED WARNINGS MAY LEAD TO ELECTRIC SHOCK, DAMAGES/LIQUID LEAKAGE, OVERHEAT, EXPLOSIONS OR FIRE AND COULD PROVOKE SERIOUS INJURIES TO PEOPLE AND/OR OBJECTS.

TREATMENT OF END-OF-LIFE BATTERIES / ACCUMULATORS

BATTERIES/ACCUMULATORS AND THE ELECTRIC OR ELECTRONIC DEVICE CANNOT BE DISPOSED AS DOMESTIC WASTE, THIS HAS TO BE DONE IN CONFORMITY WITH 2002/96/EC AND 2006/66/EC DIRECTIVES AND TO NATIONAL LAWS ON WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT. SEPARATE COLLECTION IS ESSENTIAL TO AVOID DAMAGES TO THE ENVIRONMENT.

IN ORDER TO REMOVE THE BATTERY, PLEASE REFER TO THE SPECIFIC INSTRUCTIONS.

FAILURE TO RESPECT WHAT STATED ABOVE IS SANCTION BY THE LAW.

REPLACING BATTERIES

Step 1. Open the battery compartment on the unit and remove exhausted batteries.

Step 2. Insert the batteries, making sure that the positive and negative poles are orientated according to the correct polarity.

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