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GRAZIE!

DISPOSAL OF THE ELECTRIC OR ELECTRONIC DEVICE AT THE END OF ITS LIFE

**INFORMATION FOR USERS ABOUT THE CORRECT DISPOSAL
OF THE PRODUCT ACCORDING TO THE EUROPEAN DIRECTIVE 2012/19/EU
AND TO THE LEGISLATIVE DECREE N. 49/2014.**

DEAR CUSTOMER,
IT IS IMPORTANT TO REMIND YOU THAT, AS A GOOD CITIZEN, YOU CANNOT DISPOSE OF THIS PRODUCT AS MUNICIPAL WASTE; YOU HAVE TO DO A SEPARATE COLLECTION. THE SYMBOL LABELLED ON THE APPLIANCE INDICATES THAT THE RUBBISH IS SUBJECT TO "SEPARATE COLLECTION". ELECTRICAL AND ELECTRONIC APPLIANCE CANNOT BE CONSIDERED AS DOMESTIC WASTE. THEREFORE, AT THE END OF ITS USEFUL LIFE, THIS PRODUCT SHOULD NOT BE DISPOSED WITH MUNICIPAL WASTES. THE USER MUST THEREFORE ASSIGN OR (HAVE COLLECTED) THE RUBBISH TO A TREATMENT FACILITY UNDER LOCAL ADMINISTRATION, OR HAND IT OVER TO THE RESELLER IN EXCHANGE FOR AN EQUIVALENT NEW PRODUCT.

THE SEPARATE COLLECTION OF THE RUBBISH AND THE SUCCESSIVE TREATMENT, RECYCLING AND DISPOSAL OPERATIONS PROMOTE PRODUCTION OF APPLIANCES MADE WITH RECYCLED MATERIALS AND REDUCE NEGATIVE EFFECTS ON HEALTH AND THE ENVIRONMENT CAUSED BY IMPROPER TREATMENT OF RUBBISH.

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- FOR A CORRECT AND SAFE USE OF THE APPLIANCE

ONLY FOR HOME USE. ATTENTION! READ THE INSTRUCTIONS GIVEN IN THIS HANDBOOK VERY CAREFULLY.

IT SHOULD BE PUT ASIDE FOR FURTHER CONSULTATION DURING THE ENTIRE LIFECYCLE OF THE PRODUCT SINCE IT GIVES IMPORTANT INDICATIONS ON THE PRODUCT'S FUNCTIONS AND SAFETY DURING INSTALLATION, USE AND MAINTENANCE.

GENERAL SAFETY INSTRUCTIONS

DANGER FOR CHILDREN!

- A POWER DEVICE IS NOT A TOY. AFTER USE, STORE IT OUT OF REACH OF CHILDREN. KEEP PACKING MATERIAL (PLASTIC BAGS, CARTON, ETC.) OUT OF REACH OF CHILDREN.
- THIS APPLIANCE IS INTENDED TO BE USED BY CHILDREN 8 YEARS OF AGE OR OLDER AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, SUPERVISED OR INSTRUCTED
- BY ADULTS RESPONSIBLE FOR THEIR SAFETY. THIS APPLIANCE IS NOT INTENDED FOR USE BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, UNLESS THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY AND HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION
- CONCERNING USE OF THE APPLIANCE BY A PERSON
- RESPONSIBLE FOR THEIR SAFETY. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION.
- CHILDREN SHOULD BE SUPERVISED TO MAKE SURE THEY DO NOT PLAY WITH THE APPLIANCE.

- DO NOT LEAVE THE APPLIANCE EXPOSED TO ATMOSPHERIC AGENTS (RAIN, FROST, SUN, ETC.).
- NEVER USE THE APPLIANCE:
 - WITH WET OR MOIST HANDS;
- THE APPLIANCE IS DESIGNED FOR DOMESTIC USE. PLACE THE APPLIANCE ONLY ON HORIZONTAL AND STABLE SURFACES, AWAY FROM SHOWERS, BATHTUB, SINKS, ETC.
- IF THE APPLIANCE IS WET OF GETS IN CONTACT WITH WATER, DO NOT USE IT AND HAVE IT CHECKED BY A SPECIALISED LABORATORY.

BATTERIES WARNINGS

- REMOVE THE BATTERIES FROM THE UNIT IF THEY ARE NOT USED FOR LONG PERIODS
- DO NOT RECHARGE NON-RECHARGEABLE BATTERIES
- RECHARGEABLE BATTERIES MUST BE REMOVED FROM THE PRODUCT BEFORE CHARGING THEM WITH AN EXTERNAL CHARGER.
- CHARGING MUST BE SUPERVISED BY AN ADULT
- INSERT BATTERIES ACCORDING TO THE CORRECT POLARITY
- INSTALLATION AND REPLACEMENT OF BATTERIES MUST BE CARRIED OUT BY AN ADULT
- DO NOT IMMERSE THE ITEM OR BATTERIES IN WATER
- DO NOT EXPOSE BATTERIES TO HEAT OR FIRE
- IF LIQUID SPILLED FROM BATTERIES COMES INTO CONTACT WITH CLOTHING OR SKIN, IMMEDIATELY WASH THE PART WITH CLEAN WATER AND SEEK MEDICAL ADVICE
- STOP USING BATTERIES IMMEDIATELY IF YOU NOTICE ANY ANOMALIES

- AVOID CONTACT OF BATTERIES WITH METAL OBJECTS, SUCH AS KEYS OR COINS
- BATTERIES MAY EXPLODE OR RELEASE TOXIC SUBSTANCES. RISK OF FIRE OR BURNS. DO NOT OPEN, CRUSH, MODIFY, DISASSEMBLE, PIERCE, HEAT OR BURN
- A BATTERY PACK IS TO BE UNDERSTOOD AS A SINGLE UNIT WHICH MUST NOT BE TORN, OPENED OR, IN ANY CASE, TAMPERED WITH
- DO NOT PUT BATTERIES IN YOUR MOUTH. IF SWALLOWED, PROMPTLY CONTACT YOUR DOCTOR OR LOCAL POISON CONTROL CENTRE
- DO NOT CAUSE SHORT CIRCUITS, TO AVOID BURNS OR FIRE
- IF BATTERIES SWELL, DISCONTINUE USE AND DISPOSE OF DEVICES IN ACCORDANCE WITH REGULATIONS
- BE VERY CAREFUL WHEN USING BATTERIES/BATTERIES/ ACCUMULATORS IN SERIES. NEVER MIX BATTERIES WITH DIFFERENT CHARACTERISTICS IN TERMS OF BRAND, CAPACITY, VOLTAGE. DO NOT MIX PARTIALLY DISCHARGED BATTERIES WITH CHARGED BATTERIES.
- ENSURE THAT BATTERIES ARE OF THE CORRECT TYPE AND VOLTAGE FOR THE APPLIANCE.

DANGER: DAMAGE AND/OR LEAKAGE!

FAILURE TO OBSERVE THE ABOVE MENTIONED WARNINGS MAY LEAD TO ELECTRIC SHOCK, DAMAGES/LIQUID LEAKAGE, OVERHEAT, EXPLOSIONS OR FIRE AND COULD PROVOKE SERIOUS INJURIES TO PEOPLE AND/OR OBJECTS.

TREATMENT OF END-OF-LIFE BATTERIES / ACCUMULATORS

BATTERIES/ACCUMULATORS AND THE ELECTRIC OR ELECTRONIC DEVICE CANNOT BE DISPOSED AS DOMESTIC WASTE, THIS HAS TO BE DONE IN CONFORMITY WITH 2002/96/EC AND 2006/66/EC DIRECTIVES AND TO NATIONAL LAWS ON WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT. SEPARATE COLLECTION IS ESSENTIAL TO AVOID DAMAGES TO THE ENVIRONMENT.

IN ORDER TO REMOVE THE BATTERY, PLEASE REFER TO THE SPECIFIC INSTRUCTIONS.

FAILURE TO RESPECT WHAT STATED ABOVE IS SANCTION BY THE LAW.

BATTERY REPLACEMENT

STEP 1. OPEN THE BATTERY COMPARTMENT LOCATED AT THE BOTTOM OF THE DEVICE AND REMOVE ANY EXHAUSTED BATTERIES.

STEP 2. INSERT THE BATTERIES MAKING SURE THE POSITIVE AND NEGATIVE POLES ARE ORIENTED ACCORDING TO THE CORRECT POLARITY.

STEP 3. CLOSE THE BACK COVER AGAIN.

DANGER!

FAILURE TO COMPLY WITH THE ABOVE-MENTIONED WARNINGS COULD CAUSE ELECTRIC SHOCK RISKS, BREAKAGES/LIQUID SPILLS, OVERHEATING, EXPLOSION OR FIRE, WHICH MAY LEAD TO SERIOUS PERSONAL INJURIES AND/OR DAMAGE TO OBJECTS.

Thanks for choosing smart body fat scale. This fully body composition scale can help you conveniently measure body weight, BMI, body fat rate, body water, , muscle mass, bone mass, BMR

IMPORTANT: Please read these instructions in this manual carefully before using the device. Failure to comply with the instructions given in the manual or using this scale in ways other than the ones mentioned in this manual may result in injury to yourself or damage to the scale.

RECORDED PARAMETERS

- Weight
- BMI
- Body fat rate
- Body water
- Muscle mass
- Bone mass
- BMR

INDICATION SYMBOLS

Lo	Low battery – please change the batteries
Err	Overload
C	An error has occurred during measurement. Step off the scales and step on them again to repeat the measuring process

! WARNING

- DO NOT use with medical implants such as pacemakers.
- DO NOT stand on the edge of the scale or jump on it.
- DO NOT overload the scale (MAX. 180kg / 396lb).
- DO NOT drop scale or droop and objects on it as this may damage the sensors.
- DO NOT immerse the scale in water or use chemical cleaning agents. Clean the scale with slightly damp cloth.
- NOT recommended for pregnant women.
- NOT recommended for infants, toddlers, and children under 10 years of age.
- Any measurement obtained using this device is for reference only and should not be considered as a medical opinion.
- Consult your primary care doctor or physician before making changes to your diet, exercise plan or physical activities.
- Always place the scale on a hard, dry and flat surface before measurement Make sure your feet are dry before stepping on the scale.
- Connect the batteries in accordance with the correct polarities indicated.
- Remove batteries if the scale is not be used for a prolonged period of time.
- The user under 6 years old, only weight function, no BMI and body fat, etc. data.
- • The user under 10 years old, only weight & BMI function, no body fat, etc. data.

CARE AND MAINTENANCE OF THE SCALE

- Always Keep the Scale in Horizontal Position.
- Please Keep the Scale Away from Damp Conditions and High Temperatures.
- Place the Scale on A Stable Surface.
- Please Make Sure that the Surface of the Scale Is Dry at All Times.

- Do Not Step on The Scale with Wet Feet to Avoid Slipping or Injuries that May Occur.
- Do not Put Heavy Objects on The Scale When Not in Use, Otherwise Damage May Occur.
- Treat the Scale with Care as It Is A Highly Accurate Electronic Instrument.
- Clean the Surface of The Scale with Clean, Damp Cloth. Please Do Not Use Chemicals or Harsh Abrasives for Cleaning Purposes.
- This Product Is not for Trade Purposes.

BATTERY ADVICE

- Keep the Batteries out of reach of Children.
- Do not mix Old & New Batteries with different compositions or of different Brands to prevent possible Battery Leakage or Explosions.
- Do not heat or deform the Batteries or expose them to fire. Accidents may occur.
- Waste batteries should not be disposed of with household wastes.
- Please consult your Local Authorities for further information about Battery Recycling.

KIND NOTICE

- Do not place the scale on carpets or similar soft surfaces to avoid compromise in Accuracy of the scale.
- Please keep still when you're weighing yourself on the scale to ensure accuracy in the readings.

SPECIFICATION

LCD display

Weight unit: kg/lb (switch the units in the bottom of the scale body)

Weight Limit: 6kg – 180kg (13.2lb – 396lb)

Account number: No limitation (each account can add 24 persons)

Baby mode: Allow you to scale the weight of baby

SUPPORT DEVICES

Apple 8.0 or Higher

Android™ 5.0 or Higher

INSTALL FITDAYS

1. Search "Fitdays" from iTunes Store or Google Play or scan the QR code below:

iTunes

DISPONIBILE SU
Google Play

2. Download and install the App on your device.

3. The App icon will appear on your phone or tablet after the installation is completed.

HOW TO USE APP

1. Install the batteries from the package or pull out the insulator if batteries have been inside already. Place the scale on the hard and flat ground.

Flat and hard surface




Carpet or uneven surface


2. To connect the device to the Smartphone turn on the wireless connection, find the scale name in the search list of devices to pair and pair it.

3. Register your own ID by email or social media account

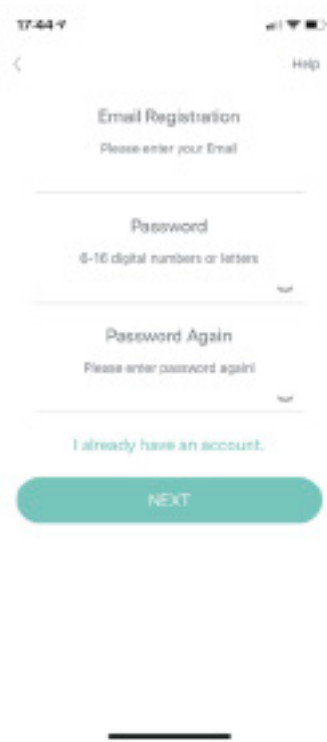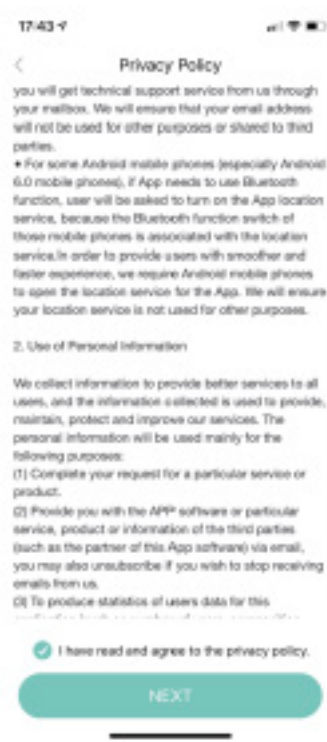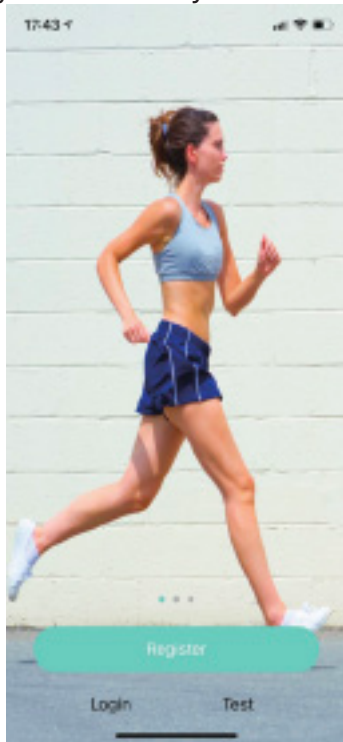4. Add personal data and confirm 0

Note: Enter the wrong gender, age, height, or corporeity level may result in measurements being inaccurate.

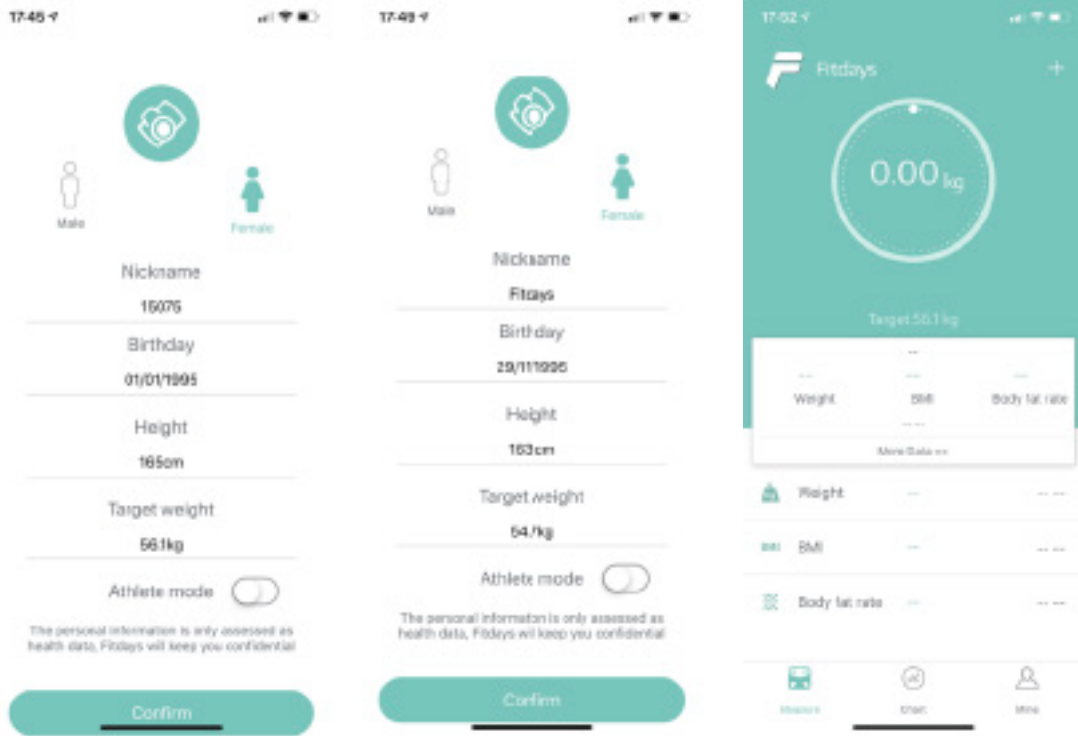5. Tap your foot on the scale to activate the display. Wait until the display shows "0.0"(kg/lb)
6. Step on the scale with BARE feet to turn the scale on.

Clean and Bare feet


With socks or shoes

NOTE: If stepping on the scale with socks, only weight and BMI will be measured and displayed on the main page. To get complete data, please step on the scale with bare feet.


7. The number on the display is flashing when the scale is weighing and analyzing. Continue to stand on the scale until the circle is stable and shows, your phone has received all the data from the scale.

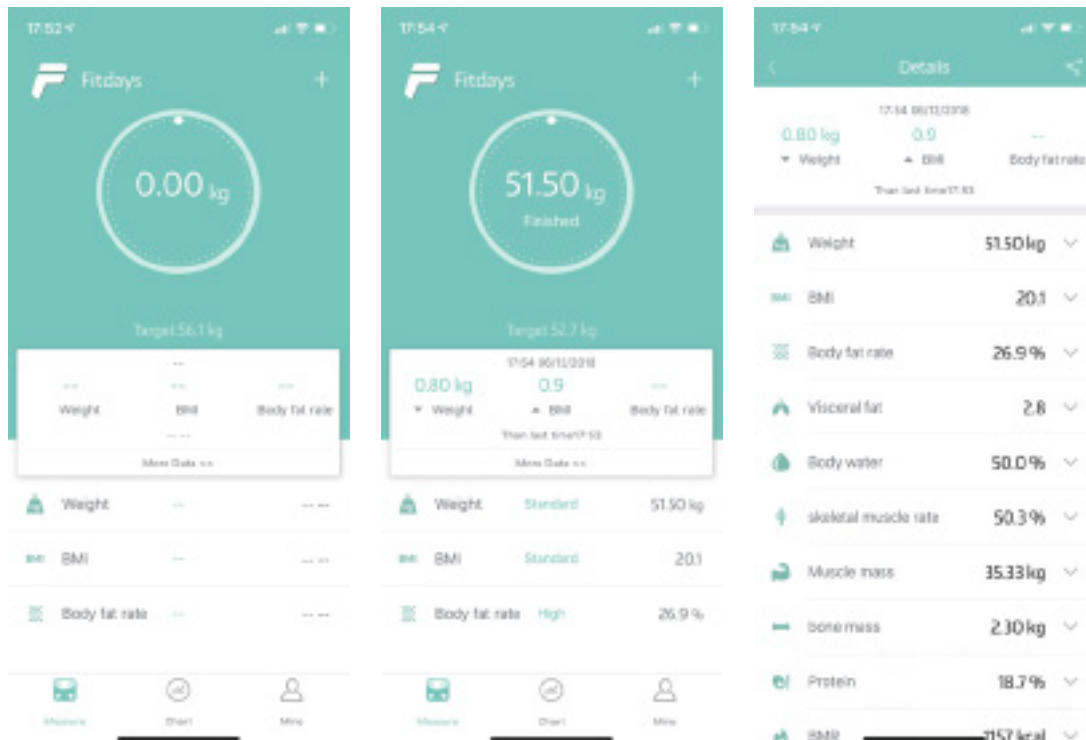8. You can find the analyzing report on Body Weight, BMI, BMR, Muscle, Bone Mass etc. (see full details on More data >>). When you click Record, you can see your long-term track chart form. Fitdays supports you to Zoom + the record for every functions via the right icon; Fitdays also supports record for recently, week, month and year;

9. Add, switch user by click  or click Me --> Manage users, then add personal data. The Smart scale can recognize which user is step on it and remind user to switch the correct user. Each account can support 24 users.

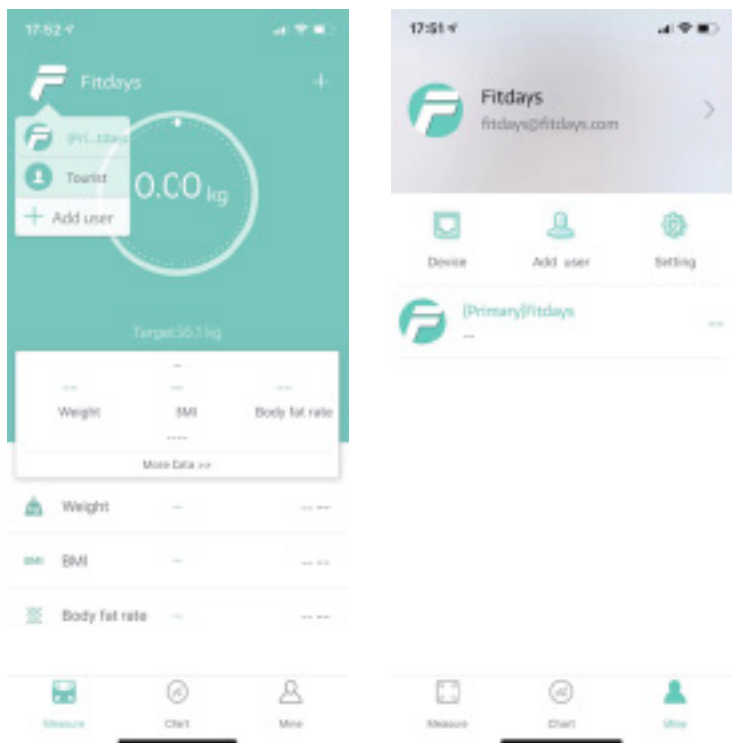10. Weighing baby: click + --> Baby Mode. Step on the scale by your own, after scale gets your weight, and then hold the baby to get the weight of baby.

TROUBLESHOOTING GUIDE

1. No Wireless connection

- The scale and smart phone must both have wireless connection on.
- Check the Apple version of your smart phone; it must be Apple 8.0 or higher / Android 5.0 or higher.
- Click Me -- > device find the scale

2. No body fat data measured when weighing

- User must be barefoot to measure the body fat.
- Make sure your feet and the scale are both dry.
- Scale not finishing weighing when the number display is flashing (Keep standing on the scale around 10 seconds until the number on the display stops flashing).
- The wireless connection is turned off.
- The device is already paired with another appliance.

3. The scale is producing inconsistent weights.

- Check whether the scale is on a hard, flat surface
- Check each sensor foot on the back to make sure nothing is stuck to the bottom of it.
- Remove the batteries to restart the scale.

4. The App is display the incorrect flat-muscle ratio

- The body parameter setting is incorrect, check to make sure that you enter the correct gender, height, and age.

5. How do I re-calibrate the scale after moving it or change the batteries?

- Step on the scale to turn it on. Let it automatically drop to 0.0kg/lb to calibrate.

6. The scale does not turn on.

- Remove the batteries and put it back. Make sure to insert the batteries in the correct direction.

7. How to change unit?

- You need to change the scale unit by the backside button. Then APP will show the same unit with the scale.

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Android, Google Play, and the Google Play logo are trademarks of Google LLC.

Apple, iPhone, iPod touch, and iTunes is a trademark of Apple Inc., registered in the U.S. and other countries.

WARRANTY CONDITIONS

This guarantee is intended to be delivered to the consumer, meaning by consumer the one who uses the good for purposes unrelated to any business or professional activity carried out, and, therefore, for domestic / private use only. Furthermore, it does not in any way affect the consumer's rights established by Legislative Decree 206/05. This certificate must be kept together with the tax receipt or other document made fiscally mandatory.

1. Galileo SpA guarantees the appliance for a period of 2 years from the purchase date (12 months if purchased with an invoice and not with a receipt).
2. The guarantee is valid exclusively on Italian territory including the Republic of San Marino and the Vatican City, and it must be proven by a receipt or any other proof of purchase, bearing the name of the Point of Sale, the purchase date of the appliance as well as the identification details of the same (type, model).
3. The guarantee covers the free replacement or repair of the component parts of the appliance which were found to be defective at the origin due to manufacturing defects.
4. Within a year of life of the product, the defect is considered as existing at the origin, subsequently its existence must be proven.
5. Knobs, lamps, glass and plastic parts, pipes and accessories, as well as all removable parts and all aesthetic parts are not covered by the guarantee, any aesthetic defects must be found and reported immediately after purchase.
6. All materials susceptible to wear (rubber / plastic seals, belts, etc.) are not covered by the guarantee unless the defect is proven to be dated back to a manufacturing origin.
7. Galileo SpA declines all responsibility for any damages that may arise, directly or indirectly to persons and things, as a consequence of failure to comply with all the instructions from the user, and especially as concerning the ones included in the specific manual for the use provided with each appliance, and, in particular, as regards the warnings relating to installation, use and maintenance of the appliance itself.

THE WARRANTY DOES NOT BECOME EFFECTIVE IN THE FOLLOWING CASES:

- a. Damage due to transport or handling of goods
- b. Incorrect installation and assembly interventions, as well as incorrect connection or incorrect supply from the electricity supply networks, lack of connection to the earth socket, and any incorrect adjustments
- c. Incorrect, imperfect or incomplete installation due to incapacity of a proper use by the User or by the installer, if the instructions for conservation and use of the same provided in the manuals are not respected, excessive environmental climatic conditions, etc.
- d. The warranty does not apply to problems exclusively related to failure to meet the customer's subjective expectations
- e. Damage caused by atmospheric agents (lightning, earthquake, fire, flood, etc.)
- f. Repairs, modifications or tampering carried out by unsuitable personnel. For failures due to lack of regular maintenance (cleaning of condenser coils, condensate water drainage holes, drain filters and water inlet, elimination of foreign bodies, limestone, etc.)
- g. No proof of purchase
- h. Improper use of the appliance, not of the domestic type or other than as established by the manufacturer
- i. For all external circumstances that are not attributable to manufacturing defects.

**Galileo S.p.A. - C.F/ P.I. 08442821008, Capitale Sociale 500.000 I.V., Rea RM - 1094899
Strada Galli n. 27, 00019 Tivoli, Roma**

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Android, Google Play e il logo di Google Play sono marchi di Google LLC.
Il logo Apple, iPhone, iPod touch, and iTunes sono marchi registrati di proprietà della Spett. Apple Inc.