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GRAZIE!

**DISPOSAL OF THE ELECTRIC OR ELECTRONIC DEVICE AT THE END OF ITS LIFE**

**INFORMATION FOR USERS ABOUT THE CORRECT DISPOSAL  
OF THE PRODUCT ACCORDING TO THE EUROPEAN DIRECTIVE 2012/19/EU  
AND TO THE LEGISLATIVE DECREE N. 49/2014.**

DEAR CUSTOMER,

IT IS IMPORTANT TO REMIND YOU THAT, AS A GOOD CITIZEN, YOU CANNOT DISPOSE OF THIS PRODUCT AS MUNICIPAL WASTE; YOU HAVE TO DO A SEPARATE COLLECTION. THE SYMBOL LABELLED ON THE APPLIANCE INDICATES THAT THE RUBBISH IS SUBJECT TO "SEPARATE COLLECTION". ELECTRICAL AND ELECTRONIC APPLIANCE CANNOT BE CONSIDERED AS DOMESTIC WASTE. THEREFORE, AT THE END OF ITS USEFUL LIFE, THIS PRODUCT SHOULD NOT BE DISPOSED WITH MUNICIPAL WASTES. THE USER MUST THEREFORE ASSIGN OR (HAVE COLLECTED) THE RUBBISH TO A TREATMENT FACILITY UNDER LOCAL ADMINISTRATION, OR HAND IT OVER TO THE RESELLER IN EXCHANGE FOR AN EQUIVALENT NEW PRODUCT. THE SEPARATE COLLECTION OF THE RUBBISH AND THE SUCCESSIVE TREATMENT, RECYCLING AND DISPOSAL OPERATIONS PROMOTE PRODUCTION OF APPLIANCES MADE WITH RECYCLED MATERIALS AND REDUCE NEGATIVE EFFECTS ON HEALTH AND THE ENVIRONMENT CAUSED BY IMPROPER TREATMENT OF RUBBISH.

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# USER INSTRUCTIONS

EN **Kooper**<sup>®</sup>

## ■ FOR A CORRECT AND SAFE USE OF THE APPLIANCE

ONLY FOR HOME USE. ATTENTION! READ THE INSTRUCTIONS GIVEN IN THIS HANDBOOK VERY CAREFULLY. IT SHOULD BE PUT ASIDE FOR FURTHER CONSULTATION DURING THE ENTIRE LIFECYCLE OF THE PRODUCT SINCE IT GIVES IMPORTANT INDICATIONS ON THE PRODUCT'S FUNCTIONS AND SAFETY DURING INSTALLATION, USE AND MAINTENANCE.

## GENERAL SAFETY INSTRUCTIONS

### DANGER FOR CHILDREN!

- A POWER DEVICE IS NOT A TOY. AFTER USE, STORE IT OUT OF REACH OF CHILDREN. KEEP PACKING MATERIAL (PLASTIC BAGS, CARTON, ETC.) OUT OF REACH OF CHILDREN.
- THIS APPLIANCE IS INTENDED FOR USE BY CHILDREN 8 YEARS OF AGE OR OLDER AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, SUPERVISED OR INSTRUCTED BY ADULTS RESPONSIBLE FOR THEIR SAFETY. THIS APPLIANCE IS NOT INTENDED FOR USE BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, UNLESS THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY AND HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY.
- CHILDREN SHOULD BE SUPERVISED TO MAKE SURE THEY DO NOT PLAY WITH THE APPLIANCE.

### DANGER OF INJURY

- PAY UTMOST ATTENTION IN CASE THE APPLIANCE IS USED IN THE PRESENCE OF CHILDREN.
- NEVER LEAVE THE APPLIANCE UNATTENDED WHILE WORKING.
- TO AVOID INJURIES, DO NOT INTRODUCE YOUR HANDS IN THE MACHINE WHILE IT IS OPERATING AND KEEP AWAY HAIR, CLOTHES AND OTHER OBJECTS, SWITCH OFF AND UNPLUG THE APPLIANCE BEFORE CHANGING THE ACCESSORIES! DO NOT CONNECT THE PLUG TO THE POWER SOCKET IF THE NECESSARY ACCESSORIES ARE NOT INSTALLED. NEVER TOUCH THE MOVING PARTS OF THE APPLIANCE AND ALWAYS

WAIT FOR THE APPLIANCE TO STOP COMPLETELY. DISCONNECT THE PLUG FROM THE OUTLET WHEN THE APPLIANCE IS NO LONGER IN USE

- DO NOT CLEAN WHILE THE APPLIANCE IS WORKING. DO NOT USE ANY SHARP OBJECT TO REMOVE POSSIBLE FOREIGN OBJECTS. THE CLEANING MUST NOT BE CARRIED OUT BY UNSUPERVISED CHILDREN

### **RISK OF ELECTRICAL SHOCK!**

- DO NOT LEAVE THE APPLIANCE EXPOSED TO ATMOSPHERIC AGENTS (RAIN, FROST, SUN, ETC.)

### **NEVER USE THE APPLIANCE:**

- IF IT HAS FALLEN AND HAS OBVIOUS RUPTURES OR MALFUNCTIONS;
- IN GAS SATURATED ENVIRONMENTS .
- REPAIRS SHOULD ONLY BE CARRIED OUT BY A QUALIFIED TECHNICIAN .
- REPAIRING THE APPLIANCE BY YOURSELF CAN BE EXTREMELY DANGEROUS .
- AN APPLIANCE THAT IS NOT WORKING ANYMORE, MUST BE RENDERED UNUSABLE, BY CUTTING THE POWER CORROBE BEFORE THROWING IT
- CLEAN THE APPLIANCE USING A DAMP CLOTH. DO NOT USE SOLVENTS OR ABRASIVE DETERGENTS .
- THE APPLIANCE IS DESIGN ED FOR DOMESTIC USE.

### **DANGER OF EXPLOSION!**

- DON'T IMMERSE THE TOOTHBRUSH BODY IN WATER
- DON'T EXPOSE TO HEAT OR FIRE IF THE LIQUID LEAKED FROM THE BATTERIES COMES INTO CONTACT WITH CLOTHES OR THE SKIN, WASH IMMEDIATELY WITH CLEAN WATER AND CONSULT A DOCTOR IMMEDIATELY. STOP IMMEDIATELY THE USE OF THE BATTERIES IN CASE OF FAULTS. AVOID CONTACT OF BATTERIES WITH METAL OBJECTS SUCH AS KEYS OR COINS
- IF YOU ARE BEING TREATED FOR ANY ORAL DISORDER, CONSULT YOUR DENTIST BEFORE USE
- THIS TOOTHBRUSH IS FOR PERSONAL USE FOR YOUR ORAL HYGIENE AND IS NOT INTENDED TO BE USED ON MULTIPLE PATIENTS IN DENTAL PRACTICES OR RELATED INSTITUTIONS. TO PREVENT BREAKAGE OF THE BRUSH HEAD, WHICH COULD RESULT IN CHOKING HAZARD DUE TO SMALL PARTS OR POSSIBLE DAMAGE TO TEETH: .
- 1. BEFORE EACH USE MAKE SURE THAT THE HEAD IS WELL INSERTED. STOP USING

THE TOOTHBRUSH IF IT IS NO LONGER POSSIBLE TO INSERT THE BRUSH HEAD CORRECTLY. NEVER USE THE BRUSH WITHOUT A BRUSH HEAD INSERTED.

- 2. IF THE BRUSH HANDLE FALLS OFF, REPLACE THE BRUSH HEAD BEFORE THE NEXT USE, EVEN IF THERE IS NO VISIBLE DAMAGE
- 3. REPLACE THE BRUSH HEAD EVERY 3 MONTHS OR SOONER IF THE BRUSH HEAD APPEARS WORN
- 4. CLEAN THE BRUSH HEAD THOROUGHLY AFTER EACH USE

**THE PRODUCT SHOULD BE USED ONLY FOR THE USE DESCRIBED IN THE MANUAL. DO NOT USE ACCESSORIES NOT RECOMMENDED BY THE MANUFACTURER.**

### **BATTERY WARNINGS**

- CHARGING MUST BE SUPERVISED BY AN ADULT
- DO NOT IMMERSE ITEM OR BATTERIES/BATTERY PACKS/ACCUMULATORS IN WATER
- DO NOT EXPOSE BATTERIES/BATTERY PACKS/ACCUMULATORS TO HEAT OR FIRE
- IF LIQUID SPILLED FROM BATTERIES/BATTERIES/BATTERY PACKS/ACCUMULATORS COMES INTO CONTACT WITH CLOTHING OR SKIN, IMMEDIATELY WASH THE PART WITH CLEAN WATER AND SEEK MEDICAL ADVICE
- STOP USING BATTERIES/BATTERIES/BATTERY PACKS/ACCUMULATORS IMMEDIATELY IF YOU NOTICE ANY ABNORMALITIES
- AVOID CONTACT OF BATTERIES/BATTERY PACKS/ACCUMULATORS WITH METAL OBJECTS, SUCH AS KEYS OR COINS
- BATTERIES/BATTERY PACKS/ACCUMULATORS MAY EXPLODE OR RELEASE TOXIC SUBSTANCES. RISK OF FIRE OR BURNS. DO NOT OPEN, CRUSH, MODIFY, DISASSEMBLE, PUNCTURE, HEAT OR BURN;
- A BATTERY PACK IS TO BE UNDERSTOOD AS A SINGLE UNIT WHICH MUST NOT BE TORN, OPENED OR, IN ANY CASE, TAMPERED WITH;
- DO NOT PUT BATTERIES/BATTERY PACKS/RECHARGEABLE BATTERIES IN YOUR MOUTH. IF SWALLOWED, PROMPTLY CONTACT YOUR DOCTOR OR LOCAL POISON CONTROL CENTRE

- DO NOT CAUSE SHORT CIRCUITS, TO AVOID BURNS OR FIRE
- IN CASE OF SWELLING OF BATTERIES/BATTERY PACKS/ACCUMULATORS, DISCONTINUE USE AND DISPOSE OF THE DEVICES IN ACCORDANCE WITH APPLICABLE REGULATIONS

NOTE: THE PRODUCT CONTAINS A NON-REMOVABLE BATTERY THAT CANNOT BE REMOVED.

**DANGER: DAMAGE AND/OR SPILLS!**

FAILURE TO OBSERVE THE FOLLOWING WARNINGS ABOVE COULD RESULT IN RISK OF ELECTRIC SHOCK, DAMAGE/OUTFLOW OF LIQUID, OVERHEATING, EXPLOSION OR FIRE AND CAUSE SERIOUS PERSONAL INJURY AND/OR PROPERTY DAMAGE.


**TREATMENT OF BATTERIES/ACCUMULATORS AT END OF LIFE.**

THE BATTERY/ACCUMULATOR/BATTERY AND ELECTRICAL AND ELECTRONIC EQUIPMENT, IF ANY, SHOULD NOT BE DISPOSED OF AS MUNICIPAL WASTE BUT IN ACCORDANCE WITH DIRECTIVES 2002/96/CE, 2006/66/CE AND NATIONAL TRANSPOSING LAWS FOR END-OF-LIFE PRODUCT DISPOSAL AND SHOULD BE SUBJECT TO SEPARATE COLLECTION TO AVOID DAMAGE TO THE ENVIRONMENT. TO REMOVE THE BATTERY, REFER TO THE SPECIFICATIONS

## FEATURES:

### Vibrating bristles movement: dynamic cleaning action up to 45,000 strokes per minute

- 4 brushing modes
- Thorough cleaning of gums and hard-to-reach areas
- 2-minute smart timer with interval pause every 30 seconds
- Replaceable toothbrush head
- Long-lasting battery: up to 3 weeks of use with a single charge
- USB charging, easy and convenient for use at home or traveling
- Charging LED indicator: low battery/full battery/charging battery.

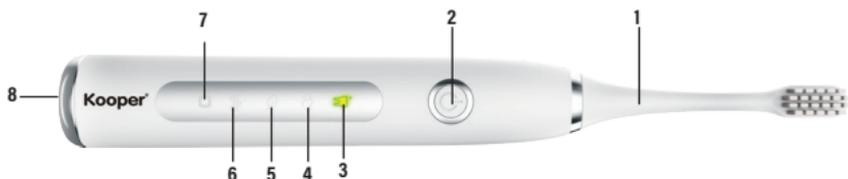1. Toothbrush head
2. On/Off button /speed/mode
3. Deep Cleaning Mode
4. Everyday cleaning mode
5. Gentle cleaning mode
6. Massage mode
7. Charging LED indicator
8. USB cable input A

### HOWTOUSE:

- Secure the toothbrush head and make sure it is locked in place
- Select the brushing mode by repeatedly pressing the power button
- Press and hold down power button 2 for 2 seconds.
- Pressing the same button sequentially switches to the various modes, keeping in mind that for each mode there are four different intensities in ascending order as shown in the following figure: (see file in folder "Electric Toothbrush Operation Modes")
- Choose the desired intensity
- Place the toothbrush on the teeth at a slight angle towards the gum line
- Press "power" button to turn it on
- Gently move the brush head over the tooth surface.
- Do not apply too much pressure on the teeth or brush an area for a long period of time. Please note it is normal to experience a vibrating sensation when starting to use an electric toothbrush. Once you get used to it, the feeling will stop
- The toothbrush will automatically turn off after 2 minutes, recommended brushing time by dentists, or by holding down the power button for 2 seconds.

**INSERT:** Slide the toothbrush head downwards to the handle, matching the grooves on the brush head **A** and the fixed part of the brush **T**.

**REMOVE:** remove gently the head **A** by sliding it in the opposite direction to the toothbrush body **T**.

**CHARGING LED LIGHTS INDICATORS:**

- **RED** charging battery/battery run out
- **GREEN** battery completely charged

**CHARGING TOOTHBRUSH:** When battery charge level is low, charging red LED indicator will light up. When battery will run out completely, **the red LED light will start to flash** and the toothbrush will turn off soon after. Plug your toothbrush by the USB cable or AC adaptor to any power source to start charging it: the indicator light will stay on. When battery is fully charged, the green indicator **light will automatically turn on**.

**Please note:** use of any other accessory may cause personal injuries or damages to things

## PRECAUTIONS :

Do not disassemble, repair or alter the product without authorization.

Persons suffering from severe periodontal disease or who have undergone oral surgery should consult a doctor before use.

Small children, infants and people who are unable to operate the device should be supervised by a parent or adult

Do not immerse the product in water

Do not leave the toothbrush in places with high temperatures or exposed to sunlight

Do not force the bristles hard against teeth or gums.

## WARRANTY CONDITIONS

This guarantee is intended to be delivered to the consumer, meaning by consumer the one who uses the good for purposes unrelated to any business or professional activity carried out, and, therefore, for domestic / private use only. Furthermore, it does not in any way affect the consumer's rights established by Legislative Decree 206/05. This certificate must be kept together with the tax receipt or other document made fiscally mandatory.

1. Galileo SpA guarantees the appliance for a period of 2 years from the purchase date (12 months if purchased with an invoice and not with a receipt).
2. The guarantee is valid exclusively on Italian territory including the Republic of San Marino and the Vatican City, and it must be proven by a receipt or any other proof of purchase, bearing the name of the Point of Sale, the purchase date of the appliance as well as the identification details of the same (type, model).
3. The guarantee covers the free replacement or repair of the component parts of the appliance which were found to be defective at the origin due to manufacturing defects.
4. Within a year of life of the product, the defect is considered as existing at the origin, subsequently its existence must be proven.
5. Knobs, lamps, glass and plastic parts, pipes and accessories, as well as all removable parts and all aesthetic parts are not covered by the guarantee, any aesthetic defects must be found and reported immediately after purchase.
6. All materials susceptible to wear (rubber / plastic seals, belts, etc.) are not covered by the guarantee unless the defect is proven to be dated back to a manufacturing origin.
7. Galileo SpA declines all responsibility for any damages that may arise, directly or indirectly to persons and things, as a consequence of failure to comply with all the instructions from the user, and especially as concerning the ones included in the specific manual for the use provided with each appliance, and, in particular, as regards the warnings relating to installation, use and maintenance of the appliance itself.

### **THE WARRANTY DOES NOT BECOME EFFECTIVE IN THE FOLLOWING CASES:**

- a. Damage due to transport or handling of goods
- b. Incorrect installation and assembly interventions, as well as incorrect connection or incorrect supply from the electricity supply networks, lack of connection to the earth socket, and any incorrect adjustments
- c. Incorrect, imperfect or incomplete installation due to incapacity of a proper use by the User or by the installer, if the instructions for conservation and use of the same provided in the manuals are not respected, excessive environmental climatic conditions, etc.
- d. The warranty does not apply to problems exclusively related to failure to meet the customer's subjective expectations
- e. Damage caused by atmospheric agents (lightning, earthquake, fire, flood, etc.)
- f. Repairs, modifications or tampering carried out by unsuitable personnel. For failures due to lack of regular maintenance (cleaning of condenser coils, condensate water drainage holes, drain filters and water inlet, elimination of foreign bodies, limestone, etc.)
- g. No proof of purchase
- h. Improper use of the appliance, not of the domestic type or other than as established by the manufacturer
- i. For all external circumstances that are not attributable to manufacturing defects.

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