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GRAZIE!

DISPOSAL OF THE ELECTRIC OR ELECTRONIC DEVICE AT THE END OF ITS LIFE

**INFORMATION FOR USERS ABOUT THE CORRECT DISPOSAL
OF THE PRODUCT ACCORDING TO THE EUROPEAN DIRECTIVE 2012/19/EU
AND TO THE LEGISLATIVE DECREE N. 49/2014.**

DEAR CUSTOMER,

IT IS IMPORTANT TO REMIND YOU THAT, AS A GOOD CITIZEN, YOU CANNOT DISPOSE OF THIS PRODUCT AS MUNICIPAL WASTE; YOU HAVE TO DO A SEPARATE COLLECTION. THE SYMBOL LABELLED ON THE APPLIANCE INDICATES THAT THE RUBBISH IS SUBJECT TO "SEPARATE COLLECTION". ELECTRICAL AND ELECTRONIC APPLIANCE CANNOT BE CONSIDERED AS DOMESTIC WASTE. THEREFORE, AT THE END OF ITS USEFUL LIFE, THIS PRODUCT SHOULD NOT BE DISPOSED WITH MUNICIPAL WASTES. THE USER MUST THEREFORE ASSIGN OR (HAVE COLLECTED) THE RUBBISH TO A TREATMENT FACILITY UNDER LOCAL ADMINISTRATION, OR HAND IT OVER TO THE RESELLER IN EXCHANGE FOR AN EQUIVALENT NEW PRODUCT.

THE SEPARATE COLLECTION OF THE RUBBISH AND THE SUCCESSIVE TREATMENT, RECYCLING AND DISPOSAL OPERATIONS PROMOTE PRODUCTION OF APPLIANCES MADE WITH RECYCLED MATERIALS AND REDUCE NEGATIVE EFFECTS ON HEALTH AND THE ENVIRONMENT CAUSED BY IMPROPER TREATMENT OF RUBBISH.

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- FOR A CORRECT AND SAFE USE OF THE APPLIANCE

ONLY FOR HOME USE. ATTENTION! READ THE INSTRUCTIONS GIVEN IN THIS HANDBOOK VERY CAREFULLY. IT SHOULD BE PUT ASIDE FOR FURTHER CONSULTATION DURING THE ENTIRE LIFECYCLE OF THE PRODUCT SINCE IT GIVES IMPORTANT INDICATIONS ON THE PRODUCT'S FUNCTIONS AND SAFETY DURING INSTALLATION, USE AND MAINTENANCE.

GENERAL SAFETY INSTRUCTIONS

DANGER FOR CHILDREN!

- A POWER DEVICE IS NOT A TOY. AFTER USE, STORE IT OUT OF REACH OF CHILDREN. KEEP PACKING MATERIAL (PLASTIC BAGS, CARTON, ETC.) OUT OF REACH OF CHILDREN.
- THIS APPLIANCE IS INTENDED FOR USE BY CHILDREN 8 YEARS OF AGE OR OLDER AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, SUPERVISED OR INSTRUCTED BY ADULTS RESPONSIBLE FOR THEIR SAFETY. THIS APPLIANCE IS NOT INTENDED FOR USE BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, UNLESS THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY AND HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY.
- CHILDREN SHOULD BE SUPERVISED TO MAKE SURE THEY DO NOT PLAY WITH THE APPLIANCE.

RISK OF INJURY!

CAUTION! FAILURE TO COMPLY WITH THE WARNINGS IN THIS MANUAL AND AN IMPROPER USE OF THE APPLIANCE MAY RESULT IN INJURY!

- BEFORE REPLACING THE ACCESSORIES, SWITCH OFF THE APPLIANCE.
- PAY UTMOST ATTENTION IN CASE THE APPLIANCE IS USED IN THE PRESENCE OF CHILDREN.
- NEVER LEAVE THE APPLIANCE UNATTENDED WHILE WORKING.
- THE APPLIANCE IS DESIGNED FOR DOMESTIC USE ONLY
- DO NOT USE THIS APPLIANCE FOR OTHER THAN ITS INTENDED ONE.
- TO AVOID INJURY AND DAMAGE, DO NOT TAMPER WITH THE DEVICE IN ANY WAY.

DANGER OF ELECTRIC SHOCK!


- DO NOT LEAVE THE APPLIANCE EXPOSED TO THE WEATHER (RAIN, FROST, SUN, ETC.).
- NEVER USE THE APPLIANCE
- WITH WET OR DAMP HANDS;
- IF IT HAS BEEN DROPPED AND HAS OBVIOUS CRACKS OR MALFUNCTIONS;
- IF THE PLUG OR POWER SUPPLY CABLE IS DEFECTIVE;
- IN DAMP OR GAS-SATURATED ENVIRONMENTS.
- ONLY PLACE THE APPLIANCE ON HORIZONTAL, STABLE SURFACES, AWAY FROM SHOWERS, BATHTUBS, SINKS, ETC.
- IF, HOWEVER, THE APPLIANCE SHOULD FALL INTO WATER, DO NOT USE IT ANY MORE, BUT HAVE IT CHECKED BY A SPECIALISED LABORATORY
- ANY REPAIRS MAY ONLY BE CARRIED OUT BY A SPECIALIST TECHNICIAN.
- REPAIRING THE APPLIANCE YOURSELF CAN BE EXTREMELY DANGEROUS.
- THE APPLIANCE IS DESIGNED FOR DOMESTIC USE. THE APPLIANCE MAY ONLY BE USED INDOOR

- DURING TREATMENT, CHECK RESULTS REGULARLY. THIS APPLIES IN PARTICULAR TO DIABETIC PEOPLE AS THEY ARE LESS SENSITIVE TO PAIN, THEREFORE RUNNING MAJOR RISKS OF INJURY.
- IN CASE OF PROLONGED USE, E.G. FOR THE REMOVAL OF CALLUSES ON THE FEET, THE DEVICE MAY HEAT UP MORE THAN NORMAL. TO AVOID SKIN BURNS, TAKE EXTENDED BREAKS BETWEEN INDIVIDUAL APPLICATIONS. FOR YOUR OWN SAFETY, CONSTANTLY CHECK THE TEMPERATURE OF THE APPLIANCE. THIS APPLIES IN PARTICULAR TO PEOPLE WHO ARE PARTICULARLY SENSITIVE TO HEAT.
- THE MANUFACTURER IS NOT LIABLE FOR DAMAGE CAUSED BY IMPROPER OR NON-COMPLIANT USE.

CLEANING AND MAINTENANCE INSTRUCTIONS

- CLEAN THE UNIT USING A DAMP CLOTH. DO NOT USE SOLVENTS OR ABRASIVE CLEANERS.
- - DO NOT MAKE ANY CLEANING ATTEMPTS WHILE THE APPLIANCE IS IN OPERATION.
- - CLEANING MUST NOT BE CARRIED OUT BY UNSUPERVISED CHILDREN.

PARTS NAME:

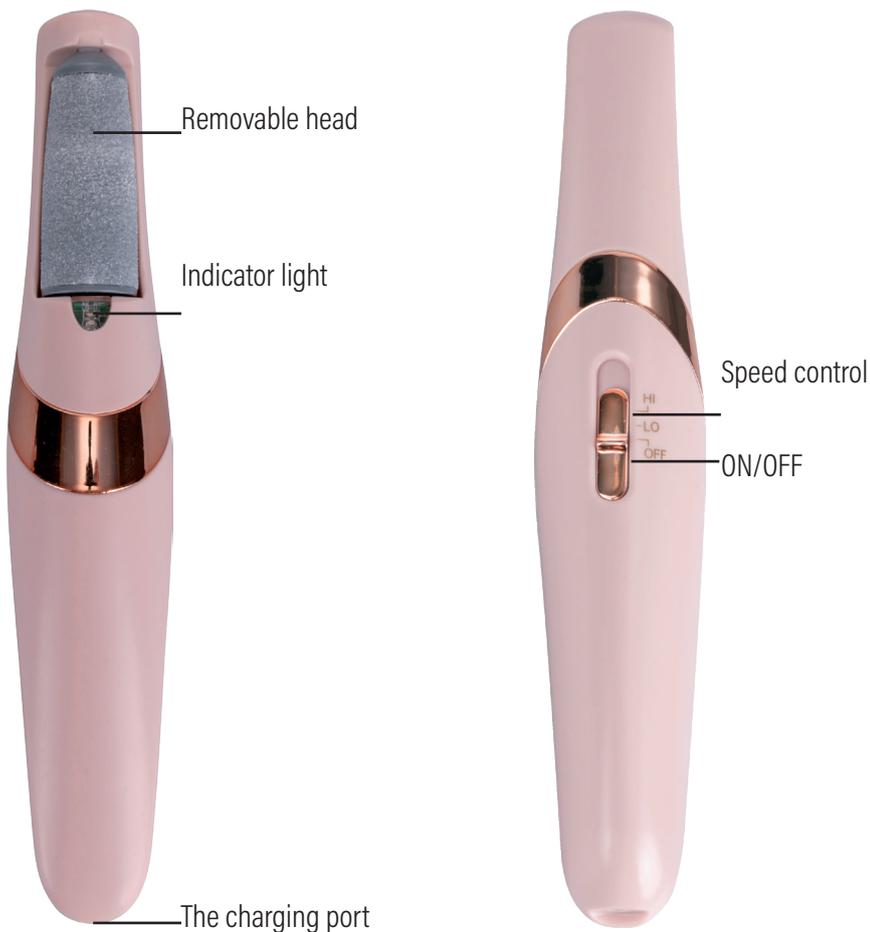INSTRUCTION FOR USE

1. Press and hold the head release button at the top of the unit and gently remove the roller head. In this way it is possible to remove and replace the desired roller head.
2. Flip the top switch upwards to turn the device on; the device operates at two speeds, select the one most suitable and least aggressive for your skin.
3. Before using the appliance, soak the feet in warm/hot water for about 20 minutes to soften the skin.
4. Switch the device on at the ON/OFF switch.
5. Gently approach the smoothing roller to the interested area without applying excessive pressure.

Nota: It is recommended to use the coarse-roller head first and then finish and smooth with the fine- roller head.

6. Wash your feet, dry them and moisturise the skin with an emollient cream..

Nota: Do not use the appliance while it is charging.

Avvertenza: the product is for domestic and non-professional use only.

PRODUCT SPECIFICATIONS:

Operating time: 20-30 minutes

Charging time: 3 hours

Battery capacity: 400 mAh

Normal temperature range for battery charging: 0-45°C

Voltage: 3.7 V

Wattage: 2 W

GENERAL WARNINGS NON-REMOVABLE BATTERIES

BATTERY WARNINGS

- Charging must be supervised by an adult
- Do not immerse item or batteries/battery packs/accumulators in water
- Do not expose batteries/battery packs/accumulators to heat or fire
- If liquid spilled from batteries/batteries/battery packs/accumulators comes into contact with clothing or skin, immediately wash the part with clean water and seek medical advice
- Stop using batteries/batteries/battery packs/accumulators immediately if you notice any abnormalities
- Avoid contact of batteries/battery packs/accumulators with metal objects, such as keys or coins
- Batteries/battery packs/accumulators may explode or release toxic substances. Risk of fire or burns. Do not open, crush, modify, disassemble, puncture, heat or burn;
- A battery pack is to be understood as a single unit which must not be torn, opened or, in any case, tampered with;
- Do not put batteries/battery packs/rechargeable batteries in your mouth. If swallowed, promptly contact your doctor or local poison control centre
- Do not cause short circuits, to avoid burns or fire
- In case of swelling of batteries/battery packs/accumulators, discontinue use and dispose of the devices in accordance with applicable regulations

Note: The product contains a non-removable battery that cannot be removed.

BATTERY WARNINGS

- Remove the battery/battery pack/accumulator from the units if they are not to be used for long periods

- Do not recharge non-rechargeable batteries
- Rechargeable batteries/battery packs/accumulators must be removed from the product before charging them with an external charger.
- Charging must be supervised by an adult
- Insert batteries/battery packs/accumulators according to the correct polarity
- Installation and replacement of batteries/battery packs/accumulators must be carried out by an adult
- Do not immerse the item or batteries/battery packs/accumulators in water
- Do not expose batteries/battery packs/accumulators to heat or fire
- If liquid spilled from batteries/battery packs/accumulators comes into contact with clothing or skin, immediately wash the part with clean water and seek medical advice
- Stop using batteries/battery packs/accumulators immediately if you notice any anomalies
- Avoid contact of batteries/battery packs/accumulators with metal objects, such as keys or coins
- Batteries/battery packs/accumulators may explode or release toxic substances. Risk of fire or burns. Do not open, crush, modify, disassemble, pierce, heat or burn
- A battery pack is to be understood as a single unit which must not be torn, opened or, in any case, tampered with
- Do not put batteries/battery packs/rechargeable batteries in your mouth. If swallowed, promptly contact your doctor or local poison control centre
- Do not cause short circuits, to avoid burns or fire
- If batteries/battery packs/accumulators swell, discontinue use and dispose of devices in accordance with regulations
- Be very careful when using batteries/batteries/accumulators in series. Never mix batteries/batteries/accumulators with different characteristics in terms of brand, capacity, voltage. Do not mix partially discharged batteries/batteries/accumulators with charged batteries/batteries/accumulators.
- Ensure that batteries/batteries/accumulators are of the correct type and voltage for the appliance.

DANGER: DAMAGE AND/OR LEAKAGE!

Failure to observe the above mentioned warnings may lead to electric shock, damages/liquid leakage, overheat, explosions or fire and could provoke serious injuries to people and/or objects.

TREATMENT OF END-OF-LIFE BATTERIES / ACCUMULATORS

Batteries/accumulators and the electric or electronic device cannot be disposed as domestic waste, this has to be done in conformity with 2002/96/ec and 2006/66/ec directives and to national laws on waste electrical and electronic equipment. Separate collection is essential to avoid damages to the environment.

In order to remove the battery, please refer to the specific instructions.

Failure to respect what stated above is sanction by the law.

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